

# A Manufacturer-Agnostic Automation Framework for SPM

Nick Sullivan Molina<sup>1</sup>, Peter Grutter<sup>2</sup>, Kirk H. Bevan<sup>1</sup>

Materials and Mining Engineering Department<sup>1</sup>, McGill University  
Physics Department<sup>2</sup>, McGill University



# Outline


1 Justification

2 Introduction to afspm


3 Design Particulars

## SPM Basics

- An **atomically-sharp tip** is **scanned** above a **surface** measuring properties.
  - Allows spatial imaging, spectroscopic analysis, sample manipulation.
- Usable on various surfaces, only require **relatively flat** sample ( $\sim 5\mu m$ ).

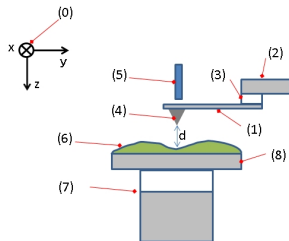

Typical AFM Configuration, Tom Toyosaki, Wikimedia Commons.

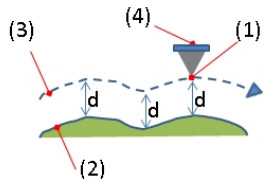

Schematics of AFM topographic image forming, Tom Toyosaki, Wikimedia Commons.

# Challenges and Automation

# Challenges and Automation

## Challenges for Wider Adoption

- 1 Preparing, running, and analyzing requires significant **domain knowledge**.
- 2 Running an experiment requires **constant user attention**.
- 3 **Statistical** understanding is **limited** by the researcher's **decisions**.

# Challenges and Automation

## Challenges for Wider Adoption

- 1 Preparing, running, and analyzing requires significant **domain knowledge**.
- 2 Running an experiment requires **constant user attention**.
- 3 **Statistical** understanding is **limited** by the researcher's **decisions**.

## Automation Prior Art

- **Conditioning** of the **tip** for proper surface characterization.
- **Structure classification**, for **where to scan** next.
- **Bayesian/Active Learning** for statistical decision making.

## A Remaining Limitation


Many SPM systems allow custom scripts to run an experiment. But:

- Scripts written for a **specific** SPM system **cannot be re-used** on others.
- While **decoupling** is possible, it is **rarely** a **priority** for researchers.

## A Remaining Limitation

Many SPM systems allow custom scripts to run an experiment. But:

- Scripts written for a **specific** SPM system **cannot be re-used** on others.
- While **decoupling** is possible, it is **rarely** a **priority** for researchers.

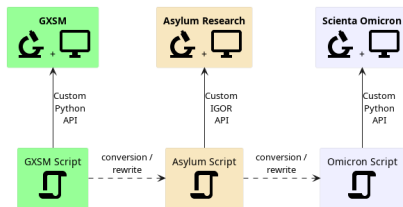




## A Remaining Limitation

Many SPM systems allow custom scripts to run an experiment. But:

- Scripts written for a **specific** SPM system **cannot be re-used** on others.
- While **decoupling** is possible, it is **rarely** a **priority** for researchers.

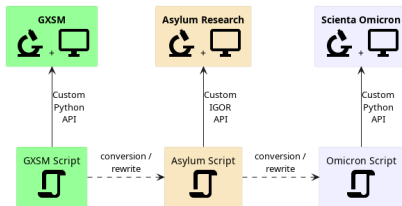

*Code **reuse** is **rare**.*

## Goals, Scope, and Design Characteristics

# Goals, Scope, and Design Characteristics

## Goals

To facilitate **code sharing** and **reusability** of developed automation.

# Goals, Scope, and Design Characteristics

## Goals

To facilitate **code sharing** and **reusability** of developed automation.

## Scope

Concerned with the **choices** a researcher would perform **during** an experiment.

# Goals, Scope, and Design Characteristics

## Goals

To facilitate **code sharing** and **reusability** of developed automation.


## Scope

Concerned with the **choices** a researcher would perform **during** an experiment.


## Design Characteristics

- **Standalone components** that each 'do one thing and do it well' (nodes).
- They communicate over '**pipes**' via generic **schemas**.
- It is easy to **split up** components among **different** computing **devices**.

# Principal Components

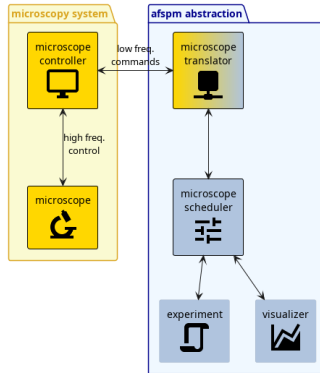


# Principal Components

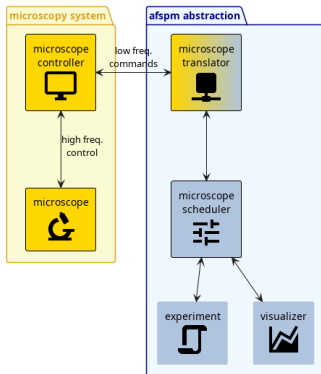

Any experiment must contain:

- **Microscope Translator:** translates between afspm-generic and microscope-specific.
- **Microscope Scheduler:** mediates control between components and translator and caches data.


# Scripting with / without afspm

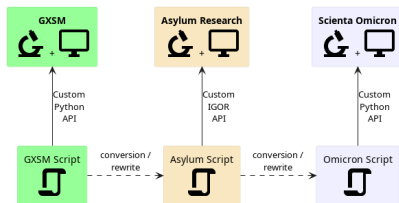


Code Reuse in SPM Scripting



# Scripting with / without afspm

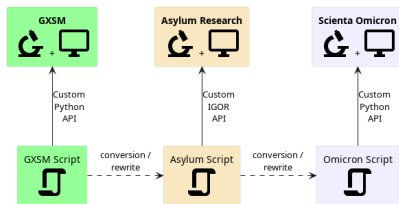


Code Reuse in SPM Scripting

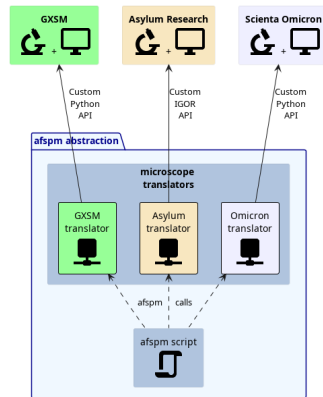


Approach Taken by afspm

# Communication Protocol

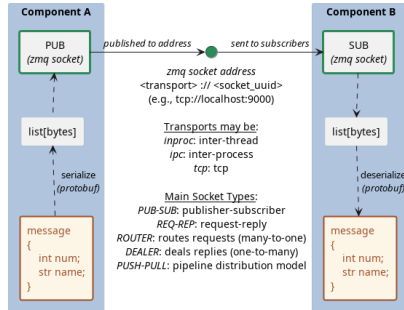

Generic *schemas* sent between *network sockets* held by components.

# Communication Protocol

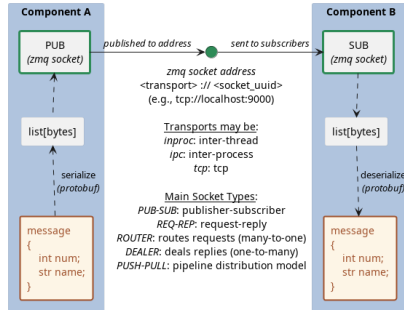


Generic *schemas* sent between *network sockets* held by components.

## Serialization / Deserialization

Translates data structures into a format that can be stored/communicated.

Google Protocol Buffers: avoids schema violations, extra error handling.

# Communication Protocol

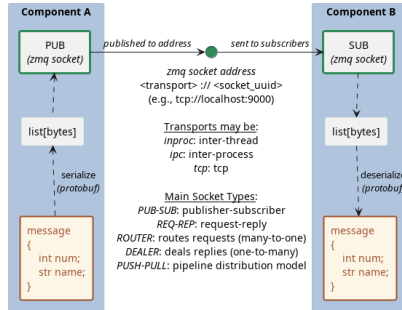

Generic *schemas* sent between *network sockets* held by components.

## Serialization / Deserialization

Translates data structures into a format that can be stored/communicated.


Google Protocol Buffers: avoids schema violations, extra error handling.

## Data Transmission Protocol


Handles sending of data between 'sockets' via common 'nodes'.

ZeroMQ: abstracts away transports used, handles common roadblocks

# Example

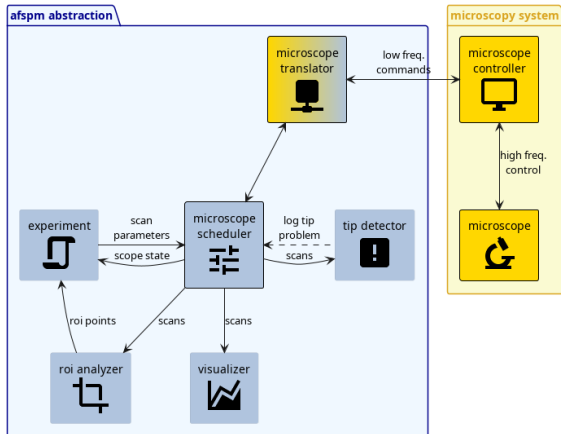

## Example

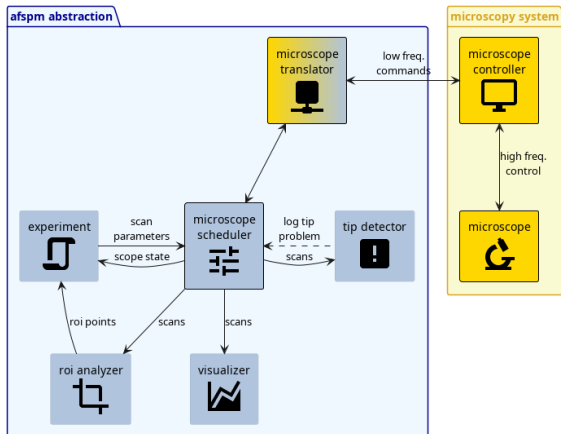

- **ROI Analyzer** reviews scans for interesting regions.
- **Experiment** switches between scanning a large region or a suggested ROI.
- **Tip Detector** evaluates the state of the tip, logs problem if deemed poor.
- **Visualizer** visualizes what is being scanned.

## Our publisher-subscriber, control-request world

We *subscribe* to events published by the SPM and *request* actions of it.

## Our publisher-subscriber, control-request world

We ***subscribe*** to events published by the SPM and ***request*** actions of it.

### Publisher-Subscriber Path

- The SPM **publishes** messages types when its 'state' changes.
- Components may **subscribe** to only messages of interest.
- Messages may be **cached** by the Scheduler, to send to new components.



# Our publisher-subscriber, control-request world

We **subscribe** to events published by the SPM and **request** actions of it.


## Publisher-Subscriber Path

- The SPM **publishes** messages types when its 'state' changes.
- Components may **subscribe** to only messages of interest.
- Messages may be **cached** by the Scheduler, to send to new components.

## Control Path

- Components send **requests** to the SPM and receive **responses**.
- The Scheduler **routes** control to the SPM (one component at a time).
- Components not in control may flag **problems** – control is dropped.
  - A user may **manually** fix a problem and unflag the problem.
  - A component that **fixes** a flagged problem can grab control.

## afspm: Detailed View

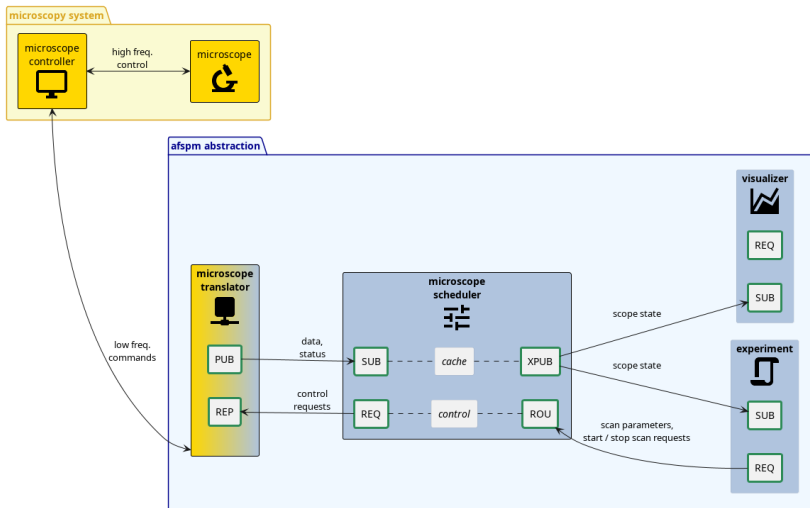

# Publisher-Subscriber Path

---

```
message DataAspects {  
  optional Size2u shape = 1;  
  optional string units = 2;  
}  
  
message SpatialAspects {  
  optional Rect2d roi = 1;  
  optional string units = 2;  
}  
  
message ScanParameters2d {  
  optional SpatialAspects spatial = 1;  
  optional DataAspects data = 2;  
}
```

---

2D Scan Parameters Schema

---

```
message Scan2d {  
  optional ScanParameters2d params = 1;  
  optional google.protobuf.Timestamp timestamp = 2;  
  optional string channel = 3;  
  optional string filename = 5;  
  repeated double values = 4;  
}
```

---

2D Scan Schema

# Publisher-Subscriber Path

---

```
message DataAspects {  
  optional Size2u shape = 1;  
  optional string units = 2;  
}  
  
message SpatialAspects {  
  optional Rect2d roi = 1;  
  optional string units = 2;  
}  
  
message ScanParameters2d {  
  optional SpatialAspects spatial = 1;  
  optional DataAspects data = 2;  
}
```

---

2D Scan Parameters Schema

---

```
message Scan2d {  
  optional ScanParameters2d params = 1;  
  optional google.protobuf.Timestamp timestamp = 2;  
  optional string channel = 3;  
  optional string filename = 5;  
  repeated double values = 4;  
}
```

---

2D Scan Schema

Frame 1

'Scan2d'

Frame 2

Scan2d Data

*Envelope*

*Serialized Data Structure*

Message Format

# Control Path

---

```
enum ControlRequest {  
    REQ_UNDEFINED = 0;  
    REQ_START_SCAN = 1;    // [...]  
    REQ_REQUEST_CTRL = 4;  // [...]  
    REQ_ADD_EXP_PRBLM = 6; // [...]  
    REQ_SET_CONTROL_MODE = 8;  
}
```

---

Control Request Schema

---

```
enum ControlResponse {  
    REP_SUCCESS = 0;  
    REP_FAILURE = 1;    // [...]  
    REP_NO_RESPONSE = 3; // [...]  
    REP_NOT_FREE = 7;    // [...]  
    REP_PARAM_NOT_SUPPORTED = 9; // [...]  
}
```

---

Control Response Schema

# Control Path

```
enum ControlRequest {  
    REQ_UNDEFINED = 0;  
    REQ_START_SCAN = 1;    // [...]  
    REQ_REQUEST_CTRL = 4;  // [...]  
    REQ_ADD_EXP_PRBLM = 6; // [...]  
    REQ_SET_CONTROL_MODE = 8;  
}
```

Control Request Schema

```
enum ExperimentProblem {  
    EP_NONE = 0;  
    EP_TIP_SHAPE_CHANGED = 1;  
    EP_DEVICE_MALFUNCTION = 2;  
    EP_FEEDBACK_NON_OPTIMAL = 3;  
}
```

Experiment Problem Schema

```
enum ControlResponse {  
    REP_SUCCESS = 0;  
    REP_FAILURE = 1;    // [...]  
    REP_NO_RESPONSE = 3;    // [...]  
    REP_NOT_FREE = 7;    // [...]  
    REP_PARAM_NOT_SUPPORTED = 9; // [...]  
}
```

Control Response Schema

```
enum ControlMode {  
    CM_UNDEFINED = 0;  
    CM_MANUAL = 1;  
    CM_AUTOMATED = 2;  
    CM_PROBLEM = 3;  
}
```

Control Mode Schema

# Control Path

```
enum ControlRequest {  
    REQ_UNDEFINED = 0;  
    REQ_START_SCAN = 1;    // [...]  
    REQ_REQUEST_CTRL = 4;   // [...]  
    REQ_ADD_EXP_PRBLM = 6;  // [...]  
    REQ_SET_CONTROL_MODE = 8;  
}
```

Control Request Schema

```
enum ExperimentProblem {  
    EP_NONE = 0;  
    EP_TIP_SHAPE_CHANGED = 1;  
    EP_DEVICE_MALFUNCTION = 2;  
    EP_FEEDBACK_NON_OPTIMAL = 3;  
}
```

Experiment Problem Schema

```
enum ControlResponse {  
    REP_SUCCESS = 0;  
    REP_FAILURE = 1;    // [...]  
    REP_NO_RESPONSE = 3;    // [...]  
    REP_NOT_FREE = 7;    // [...]  
    REP_PARAM_NOT_SUPPORTED = 9; // [...]  
}
```

Control Response Schema

```
enum ControlMode {  
    CM_UNDEFINED = 0;  
    CM_MANUAL = 1;  
    CM_AUTOMATED = 2;  
    CM_PROBLEM = 3;  
}
```

Control Mode Schema

Frame 1

REQ\_REQUEST\_CTRL

Frame 2

EP\_TIP\_SHAPE\_CHANGED

*Command Request Enum*

*Serialized Data Structure(s)*

Message Format

# Spawning, Monitoring, and the Config File

*An experiment defined via a **config** file, with components **spawned** as needed.*



# Spawning, Monitoring, and the Config File

*An experiment defined via a **config** file, with components **spawned** as needed.*

## Spawning

- Each component is instantiated as a child process (independent memory).
  - One crashed component does not crash all.
- We can spawn different components on different computers.
  - The 'spawn' command allows indication of which components to spawn.

# Spawning, Monitoring, and the Config File

*An experiment defined via a **config** file, with components **spawned** as needed.*

## Spawning

- Each component is instantiated as a child process (independent memory).
  - One crashed component does not crash all.
- We can spawn different components on different computers.
  - The 'spawn' command allows indication of which components to spawn.

## Monitoring

- The parent process monitors spawned components for heartbeats.
- If a process has not 'beat' in a timeframe, it is killed and respawned.
- Then the caching logic chosen is important – must be able to run.

# Sample Config File

---

```
# --- URLs --- #
pub_url = "tcp://127.0.0.1:9000"
psc_url = "tcp://127.0.0.1:9001"

server_url = "tcp://127.0.0.1:6666"
router_url = "tcp://127.0.0.1:6667"
# [...]
exp_scan_res = [256, 256]
# [...]
```

---

## General Variables

---

```
# ----- PubSub ----- #
[translator_pub]
class = 'afspm.io.pubsub.publisher.Publisher'
url = 'pub_url'
# [...]
[experiment_sub]
class = 'afspm.io.pubsub.subscriber.Subscriber'
sub_url = 'psc_url'
# [...]
```

---

## Intermediary Classes

---

```
[translator]
component = true
class = 'afspm.components.microscope.translators.gxsm.translator.GxsmTranslator'
publisher = 'translator_pub'
control_server = 'translator_server'
# [...]
```

---

## Config Components

## Expandability

*We defined a '**basis**' set of schemas we believe are standard.  
But this is quite arbitrary! Here are ways you can expand upon them.*

# Expandability

*We defined a '**basis**' set of schemas we believe are standard.  
But this is quite arbitrary! Here are ways you can expand upon them.*

## Operating Modes

- Defined as an SPM config describing cantilever control during a scan.
- We assume experiment has scan and spectroscopy modes configured.
- Can change by providing **op. mode** id and <key:value> **params** map.

# Expandability

*We defined a '**basis**' set of schemas we believe are standard.  
But this is quite arbitrary! Here are ways you can expand upon them.*

## Operating Modes

- Defined as an SPM config describing cantilever control during a scan.
- We assume experiment has scan and spectroscopy modes configured.
- Can change by providing **op. mode** id and <key:value> **params** map.

## Actions

- Many manufacturers provide some automation tasks (e.g. approach tip).
- Can send these by providing **action** id and <key:value> **params** map.

# Expandability

*We defined a '**basis**' set of schemas we believe are standard.  
But this is quite arbitrary! Here are ways you can expand upon them.*

## Operating Modes

- Defined as an SPM config describing cantilever control during a scan.
- We assume experiment has scan and spectroscopy modes configured.
- Can change by providing **op. mode** id and <key:value> **params** map.

## Actions

- Many manufacturers provide some automation tasks (e.g. approach tip).
- Can send these by providing **action** id and <key:value> **params** map.

## Custom Messages

- Any user may create custom messages or expand existing ones.

# The End

Let us know what you think and help us make it better.

[afspm on github](#)